

# BRSReadout Documentation

M.P.Ross, K.Venkateswara

March 12, 2019

## 1 Introduction

This document describes the “BRSReadout” software which was developed by the Eöt-Wash group as the readout software of the Beam Rotation Sensor (BRS). The BRS’s angular readout is achieved using a multi-slit autocollimator refined from the device described in [1]. The goal of the software is to turn line camera images into an number of angular readout channels, which are sent to the accompanying TwinCAT code, while also giving the user an interface to view the camera images, plot the angular data, record data to disk, and monitor performance. The source for BRSReadout software can be found here: <https://github.com/mpross/BRSReadout> or [slowcontrols/trunk/BRS2C#/BRSReadout/](https://github.com/mpross/BRS-Controls). The accompanying TwinCAT code can be found at <https://github.com/mpross/BRS-Controls> or [slowcontrols/trunk/TwinCAT3/BRS/](https://github.com/mpross/BRS-Controls) with a version for each BRS deployed at the observatories. The executable can be found under `BRSReadout/BRSReadout/bin/Debug/BRSReadout.exe`.


## 2 Installation

Once the software is downloaded either through git or svn, there are a few minor steps that need to be completed to finished the installation. For each device there is a install batch script in the main directory which replaces the configuration file (`BRSReadout/BRSReadout/bin/Debug/BRSReadout.exe.config`) with the file corresponding to the given device held in the `configs/` directory. It additionally adds a shortcut on the desktop to `BRSReadout/BRSReadout/bin/Debug/BRSReadout.exe` for convenience.

## 3 Overview


The code is split between a number of simultaneous threads to allow for increased functionality while maintaining performance. The core of the software is the `Form1.cs` class which analyzes the images as described in the Image Analysis section, calculates the various signals and status bits, and handles all of the data communication with other threads and the accompanying TwinCAT code. The graphing (`Form2.cs`), frame capturing (`Camera.cs`), and data writing (`DataWriting.cs`) are all handled in their own threads to allow the various subroutines to run at different rates while also decoupling the performance of the threads.

Upon initialization, the main loop starts the other threads, reads a collection of settings from the configuration files, and triggers the camera to begin continuous acquisition. After which is waits for a frame to come into the queue before running the frame through the analysis, calculating signals, then sending the signals to the other threads and the TwinCAT code. It then waits for another frame and repeats until the software is closed. It should be noted that the analysis is highly CPU intensive and requires a dedicated core to achieve maximal rate while avoiding interrupts from other processes. This is achieved by utilizing the CPU Affinity setting in Windows which is set using a the `AffinitySet.bat` script on startup.

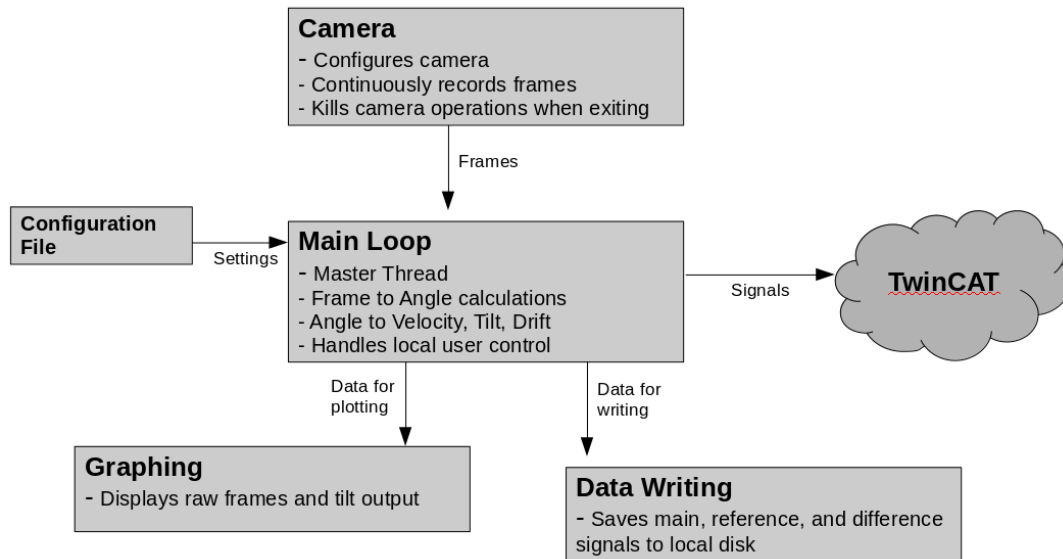

## 4 Image Analysis

The multi-slit autocollimator shines light from a fiber coupled LED through a mask consisting of 38 slits which is then set through a beam splitter. The reflected beam is imaged onto a CCD line camera while the transmitted beam reflects off of a mirror that is attached to the beam before also being imaged on the CCD. This forms two nearly identical patterns, the reference and main patterns which respectively correspond to the reflected and transmitted path of the beam splitter. The major differences between the patterns is the height of the peaks which is due to the main pattern passing through the beam splitter twice and the spatial separation of the two patterns which is due to a static angle between the main mirror and the optical path. Introducing this static angle also reflects away the secondary reflections out of the imaging plane of the CCD.

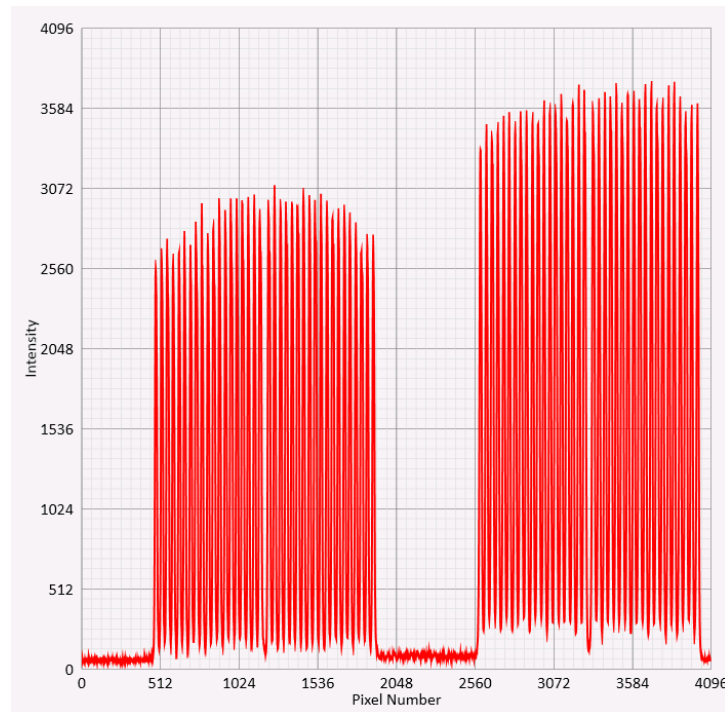

As the angle between the beam and the housing change, the main pattern (shown above on the left) translates across the CCD while the reference pattern stays stationary. The program now needs to take these images and turn it a data stream of angular readings by measuring the distance between the patterns for each frame.

The frames are recorded by the camera as a 4096 element vector of 12 bit numbers representing the intensity of each CCD pixel. The control of the camera and initial frame handling is done in the Camera.cs class which pulls the frames from

the camera using the Basler Pylon APIs, averages a small number of frames together to decrease downstream load, and puts them in a queue to be ready by the main loop. Additionally, for some of the BRSs the cameras were intentionally installed backwards. To account for this each frame is reversed before being sent on.

---

#### Frame recording code (Camera.startFrameGrab):

---

```
while (true)
{
    if (cameraType == "basler")
    {
        PylonBuffer<Byte> buffer;
        Pylon.WaitObjectWait(hWait, 100);

        Pylon.StreamGrabberRetrieveResult(hGrabber, out grabResult);

        bufferIndex = (int)grabResult.Context;

        while (buffers.TryGetValue(grabResult.hBuffer, out buffer) != true) ;
        Buffer.BlockCopy(buffer.Array, 0, data, 0, buffer.Array.Length);

        if (frameReverse)
        {
            Array.Reverse(data);
        }
        frameAveraging(data);
        Pylon.StreamGrabberQueueBuffer(hGrabber, grabResult.hBuffer, bufferIndex);
    }
}
```

---

Once a frame is placed in the queue the main loop (Form1.cs) retrieves it and under goes a series of steps to extract the position of each pattern. The first step of this process is to split the frame in half and designating one side at the reference pattern and the other as the main pattern. This splitting is done at a static location that can be set by changing the "splitPixel" parameter in the configuration file and is normally set to be as close to the reference pattern as possible. This is to give the main pattern, which will be moving during operation, the most room possible. From here on the patterns are treated as separate and undergo identical analyses to yield two channels, the main channel and the reference channel, which correspond to the location of the respective patterns.

To give the algorithm a frame to compare with, the first frame that comes into the algorithm is saved and treated as a standard against which all future frames will be compared. Since this frame is captured on the given instrument, it will contain similar nuance features as future frames and thus be a better standard than any theoretical pattern. This frame is stored and used until either the software is restarted or the "Recapture Frame" button is pressing. This recapture function is only necessary if the first frame is corrupted due to a bad camera frame or if either pattern is out of range.

Once another frame comes in, the algorithm must extract the current position of the two patterns with respective to their position in the first frame. This is done by combining two methods: a course threshold based edge finder and a fine cross correlation based center calculation. The edge finder finds the first pixel that has an intensity greater than the threshold that is set by the "threshold" parameter in the configuration file. A number set by the "pixelMargin" parameter is subtracted from this pixel number to ensure that the entire pattern is used in the cross correlation step. This is then saved as the "startIndex" and provides a measure of the patterns location to within a pixel.

---

#### Edge finding code (Form1.Pattern):

---

```
for (int j = splitPixel; j < frame.Length; j++)
{
    if (frame[j] > threshold)
    {
        startIndexRight = j - pixelMargin;
        break;
    }
    if (j == frame.Length-1)
    {
        startIndexRight = 0;
    }
}
```

---

To get sub-pixel resolution of the position, we take the cross correlation of the section of the current frame from the start index to the start index plus the pattern length (which is set by the "patternLength" parameter) and the corresponding

section in the first frame. This gives us a oscillating pattern with a central peak that is displaced from zero by the distance between the patterns. A gaussian is then fit to a number of points around this central peak. The mean value of this gaussian then corresponds to the sub-pixel displacement of the pattern and it's respective first frame pattern. This is then added to the start index to yield an accurate measure of the patterns location.

#### Cross correlation code (Form1.Pattern):

---

```
for (int k = -halflength; k <= halflength; k++)
{
    sum = 0;
    for (int m = 0; m < length; m++)
    {
        if ((m + startIndexRight + k) > 0 && (m + startIndexRightRef) > 0)
        {
            if ((m + startIndexRight + k) < frame.Length
                && (m + startIndexRightRef) < refFrame.Length)
            {
                sum += frame[m + startIndexRight + k]
                    * refFrame[m + startIndexRightRef];
            }
        }
    }
    crossCor[k + halflength] = sum;
}

//Sums x,x^2,x^3,x^4,ln(y),x ln(y),x^2 ln(y)
y = 0;
xy = 0;
xxy = 0;

for (int j = 0; j < fitLength; j++)
{
    y += Math.Log(crossCor[j + 1]);
    xy += j * Math.Log(crossCor[j + 1]);
    xxy += j * j * Math.Log(crossCor[j + 1]);
}

//Solves system of equations using Cramer's rule
D = N * (xSquar * xFourth - xCube * xCube)
    - x * (x * xFourth - xCube * xSquar) + xSquar * (x * xCube - xSquar * xSquar);
//Da = y * (xSquar * xFourth - xCube * xCube)
    - x * (xy * xFourth - xCube * xxy) + xSquar * (xy * xCube - xSquar * xxy);
Db = N * (xy * xFourth - xCube * xxy)
    - y * (x * xFourth - xCube * xSquar) + xSquar * (x * xxy - xy * xSquar);
Dc = N * (xSquar * xxy - xy * xCube)
    - x * (x * xxy - xy * xSquar) + y * (x * xCube - xSquar * xSquar);
//a = Da / D;
b = Db / D;
c = Dc / D;

mu = -b / (2 * c);


newdata[frameNo, 0] = mu + startIndexRight;
timestamps[frameNo] = data.TimeStamp(frameNo);

angleLastValue = mu + startIndexRight;
```

---

Once both patterns go through this procedure, the algorithm outputs two data streams: the raw reference pattern motion as a measure of the instruments readout noise and the difference between the two patterns which will be proportional to the tilt with decreased contribution of the readout noise due to the subtraction. These two channels are then sent on to the graphing thread, the data writing thread, and to the signal processing to be calibrated, conditioned, and split into the various output channels.

A demonstration of this analysis can be seen in the accompanying BRSReadoutDescription.m Matlab script.

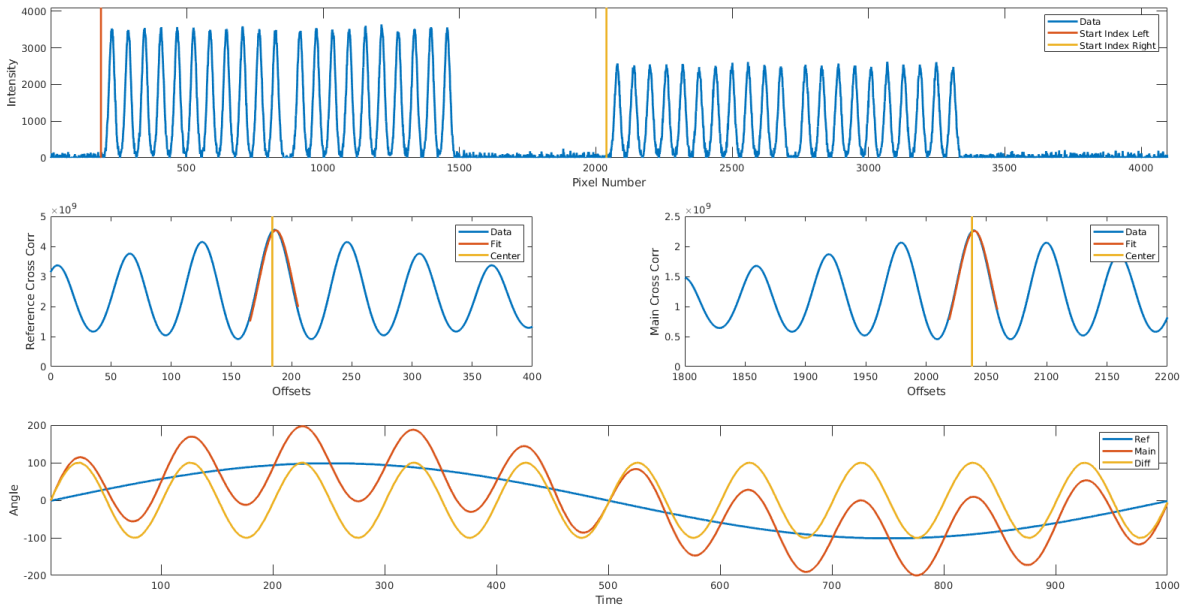

## 5 Signal Processing

Under normal conditions the raw angle (difference) signal can take on a significant range of values. Since these signals will be sent to the seismic front ends using DACs, using only the raw signal would lead to significant bit noise. To circumvent this the angle signal is split into two: a drift channel which is raw signal with a gain which matches the dynamic range of the autocollimator to the range of the DACs and a tilt signal which passed through a roughly 1 mHz high pass filter before being calibrated to 1 nrad per DAC count. The tilt channel is then a low noise signal that can be sent on to the control systems while the drift channels allow for a diagnostic monitor that ensures that the instrument is in range while being used. To allow the device to output a drift channel while the main pattern is off the edge of the CCD, the drift channel can be rescaled and recentered by triggering the "Drift override" state. **This compromises the fidelity of the tilt channel and should only be used during diagnostics.**

In addition, the instrument has two capacitor plates that are used to damp the amplitude of the beam. This requires the velocity of the angular motion around the resonant frequency of the beam which is calculated by band passing the tilt signal around the resonance (3-8 mHz) and then taking the difference. This is then square rooted to account for the fact that the force between capacitor plates goes as the voltage squared and then limited to be within the DACs range.

The tilt, drift, velocity, and reference channels are then sent out to the accompanying TwinCAT code in addition to a selection of status bits. These give the status of whether the Readout code is running, whether the camera is getting enough light, and whether the communication with the camera is functioning. Once these signals are sent the main loop cycles through for another frame until the software is shutdown.

### Signal processing code (Form1.dataSend):

---

```
angle = data[0];
refAng = refGain * (data[1] - refZeroValue);
//Has a DC subtraction to help filters. Should be about the center of the signal.
if (firstValueCounter < 40)
{
    firstValueCounter++;
    zeroValue = data[0];
    refZeroValue = data[1];
}
xHighPass[0] = angle - zeroValue;
xBandLowPass[0] = angle - zeroValue;

//Drift signal calculation, just scaled signal
if (dampOverride)
{

```

```

        drift = driftOverGain * (angle - driftOverOffset);
    }
    else
    {
        drift = driftGain * (angle - driftOffset);
    }
    //Drift rail logic
    if (Math.Abs(drift) >= 32760)
    {
        drift = Math.Sign(drift) * 32760;
    }

    //Tilt signal calculations, high pass at 10^-3 Hz then scaling.
    yHighPass[0] = highCoeff[3] * yHighPass[1] + highCoeff[4] * yHighPass[2]
        + highCoeff[0] * xHighPass[0] + highCoeff[1] * xHighPass[1]
        + highCoeff[2] * xHighPass[2];
    tilt = angleGain * yHighPass[0];

    //Capacitor signal calculations, low passed at Hz then high passed
    at Hz then differentiated and scaled
    yBandLowPass[0] = bandLowCoeff[3] * yBandLowPass[1] + bandLowCoeff[4] * yBandLowPass[2]
        + bandLowCoeff[0] * xBandLowPass[0] + bandLowCoeff[1] * xBandLowPass[1]
        + bandLowCoeff[2] * xBandLowPass[2];
    xBandHighPass[0] = yBandLowPass[0];
    yBandHighPass[0] = bandHighCoeff[3] * yBandHighPass[1] + bandHighCoeff[4] * yBandHighPass[2]
        + bandHighCoeff[0] * xBandHighPass[0] + bandHighCoeff[1] * xBandHighPass[1]
        + bandHighCoeff[2] * xBandHighPass[2];
    velocity = velGain * Math.Round((double)200 / bufferSize) *
        (yBandHighPass[0] - yBandHighPass[1]);

    if (velocity > 0)
    {
        velocity = Math.Sqrt(velocity);
    }
    else
    {
        velocity = -Math.Sqrt(-velocity);
    }


    if (Math.Abs(velocity) >= 30760)
    {
        velocity = Math.Sign(velocity) * 30760;
    }
    if (twinCatBool)
    {
        ds = new AdsStream(28);
        BinaryWriter bw = new BinaryWriter(ds);
        //Tilt signal. TwinCAT variable tilt at MW0.
        bw.Write((int)tilt);
        //Drift signal. TwinCAT variable drift at MW1
        bw.Write((int)drift);
        //Velocity signal. TwinCAT variable cap at MW2
        bw.Write((int)velocity);
        voltagerwrite = velocity / 3276;
        //Reference signal. TwinCAT variable ref at MW3
        bw.Write((int)refAng);
        //C# pulse. Sets TwinCAT variable cPulse=1 at MW4
        bw.Write((int)1);
        //Light source status bit at MW5
        bw.Write((int)lightSourceStatus);
        //Camera status bit at MW6
        bw.Write((int)cameraStatus);
        tcAds.Write(0x4020, 0, ds);
    }

```


## 6 User Interface

The user interface is design to mainly to allow a local user to monitor the status of the software and device while also allow some limited user controls. It is split between two separate windows, "BRS Readout" and "BRS Graph", corresponding to the Form1.cs and Form2.cs classes. These run in independent threads to decouple the performance of the graphing from the main loop while also utilizing multiple cores.

"BRS Readout" is the main window which will be open whenever the software is running and holds all of the user controls. It displays four statistics at the top of the panel: time since the code was run in seconds, the length of the camera queue and the graphing queue, and the voltage being written to the capacitor plates. These are used as diagnostics to ensure that the code is running as expected. The panel also has five buttons which give the user the ability to record the data to the local disk, clear the angle graph, recapture the first frame, open the graphing window, and to turn on the damping override. The ability to record to disk does not affect the output data streams but is intended to allow the BRS to collect data without communicating with the LIGO systems. The frame recapture function can be used if the first frame captured by the camera is corrupted while the damping override will allow the BRS to damp the beam motion even if the main pattern is slightly off the CCD. **This compromises the fidelity of the tilt channel and should only be used during diagnostics.**

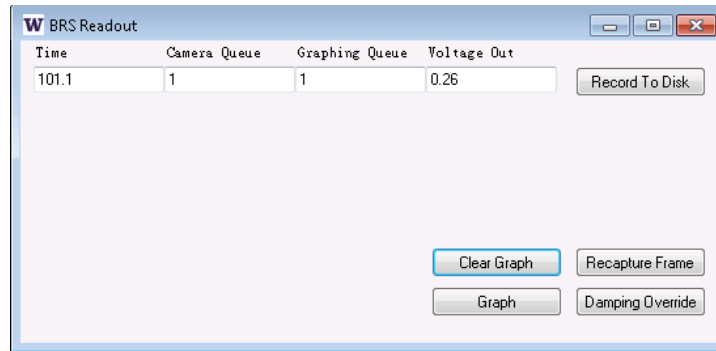

The "BRS Graph" window can be opened using the "Graph" button and holds two panes which show plots of the current frame on the left and a time series of the angular output on the right. This allows the user to get a image of the current patterns while also monitoring performance of the tilt channel. The angular plot can be cleared by using the "Clear Graph" button and plots the number of points set by the "numberOfGraphPoints" configuration parameter. To decrease the CPU load of this thread, the graphs are updated after a number of frames set by the "graphingFrameNumber" parameter. It should also be noted that the graph can experience glitches when the CPU load is too high such as when a user begins a remote desktop instance.

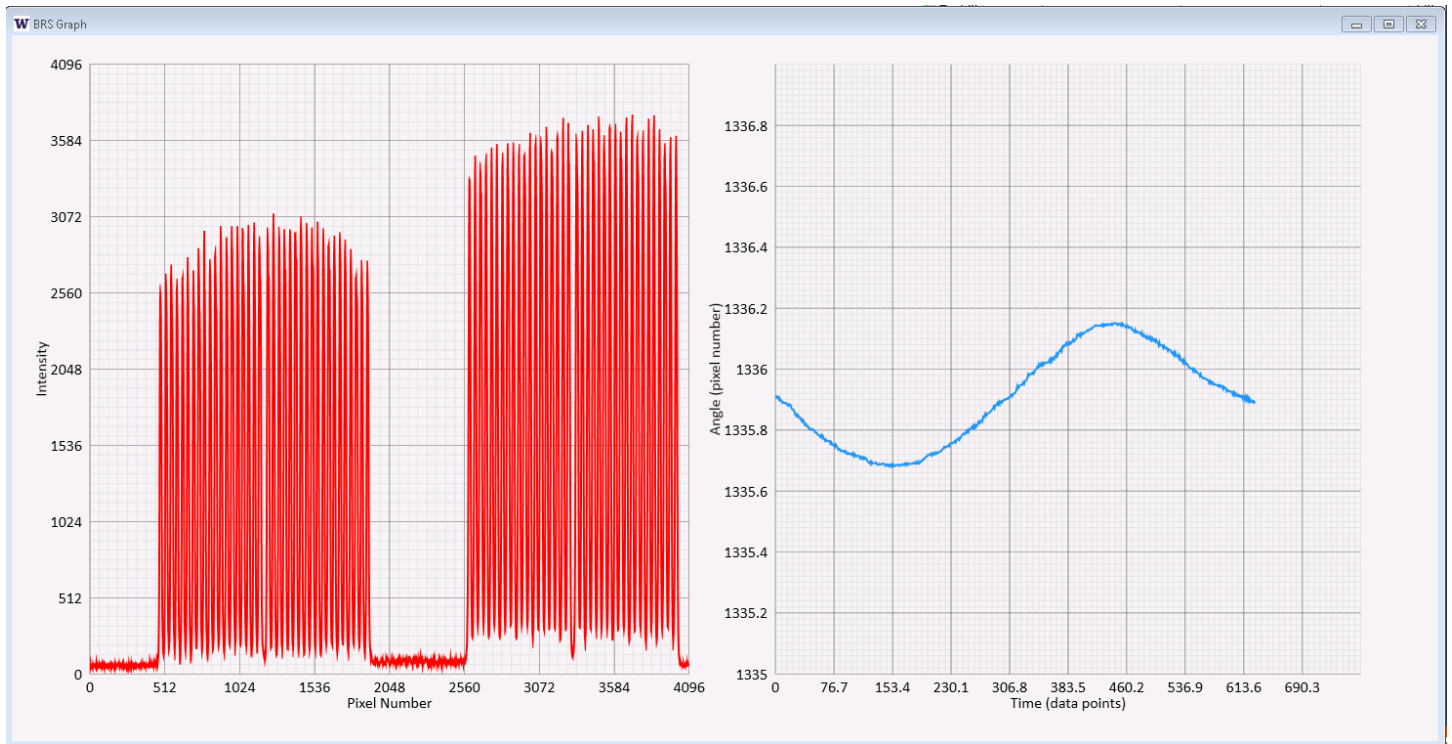

## 7 Configuration File

The parameters for each instance of the software are stored in a configuration file under `BRSReadout/BRSReadout/bin/Debug/BRSReadout.exe.config` and are loaded into the software on startup. These files can be edited manually using a standard text editor and are also initialized with default setting during the install process. The default configuration files for each device are held under `configs/`. Each parameter has a line such as `<add key="name" value="value"/>` which is imported with the following: `ConfigurationManager.AppSettings.Get("name")`. This outputs a string which then must be parsed into the correct data type for the given variable.

### General settings

- location - Device location name. Example: LHO End-X
- twinCat - Boolean controlling use of TwinCAT software. Default: true

### Camera settings

- camera - Camera type. Default: basler
- cameraBitDepth - Camera bit depth. Default: 12
- cameraExposureTime - Camera exposure time. Default: 500
- cameraWidth - Camera pixel number. Default: 4096
- ipAddress - IP Address of camera. Only needed for cameras that are routed through switch.
- frameAverageNum - Number of frames to average within Camera class. Default: 4
- frameReverse - Frame flip boolean. Only true for devices with camera installed backwards.

### Graphing settings

- graphingFrameNumber - Number of frames to wait until graph is updated. Default: 5
- bufferSize - Number on frames in buffer. Default: 20
- numberOfGraphPoints - Number of points to graph. Default: 1000

### Email settings



- emailList - List of emails to send crash reports to separated by ",".

#### Filter settings

- angleLowPass - Cutoff frequency for tilt low pass in Hz. Default: 0.1
- angleHighPass - Cutoff frequency for tilt high pass in Hz. Default: 0.0005
- velocityLowPass - Cutoff frequency for velocity low pass in Hz. Default: 0.02
- velocityHighPass - Cutoff frequency for velocity high pass in Hz. Default: 0.002

#### Fitting settings

- threshold - Edge finder threshold. Default: 1200
- splitPixel - Pixel on which frame is split between patterns. Should be set as close to the right edge of the reference pattern as possible Default: 1800
- pixelMargin - Margin to shift back by in edge finding to ensure entire pattern is captured. Default: 40
- patternLength - Pattern length. Default: 1500

#### Gain settings

- refGain - Reference channel gain. Default: 3499.8
- angleGain - Angle channel gain. Default: 3645
- velGain - Velocity channel gain. Default: -5000000000
- driftGain - Drift channel gain. Default: 60
- driftOffset - Drift channel zero pixel. Default: 2100
- driftOverGain - Drift channel gain during override. Default: 30
- driftOverOffset - Drift channel zero pixel during override. Default: 2600

## 8 Troubleshooting

The software is designed to run with very little user input however there are a number of faults which have been observed that have simple solutions.

- Angular output is outputting pure noise
  1. Ensure patterns are not overlapping by inspecting the pattern pane in the graphing window. If overlapping either recenter the beam balance or wait for it to drift back into range.
  2. Try recapturing the first frame. This will fix the issue of a corrupted first frame and should take immediate effect
- Angular output has sporadic glitches
  1. Ensure patterns are not overlapping by inspecting the pattern pane in the graphing window. If overlapping either recenter the beam balance or wait for it to drift back into range.
  2. Check the queue numbers in the main window which should both be at either 1 or 0. These will rise when the CPU is heavily loaded which will cause the program to halt momentarily causing glitches. Either decrease the CPU load by exited spurious programs or by tweaking the "graphingFrameNumber" or "frameAverageNum" parameters.
- Angular output flat lines
  1. Ensure that the main pattern is in range by inspecting the pattern pane in the graphing window. If not either recenter the beam balance or wait for it to drift back into range.
  2. Check that the time in the main window is ticking. If not the code is frozen and will likely need to be restarted.
  3. Check the "splitPixel" parameter. It should be set to a pixel that is close to the reference pattern. If set too high the working range of the instrument will be decreased. To fix decrease the parameter to roughly ten pixels above edge of the reference pattern, save configuration file, and restart code.

## 9 Previous Versions

Although this document was written specifically about the "BRSReadout" software, previous version such as the "FringeAutocollimator" software follow an almost identical image analysis and have similar GUI components.

## References

- [1] T B Arp, C A Hagedorn, Stephan Schlamming, and J H Gundlach. A reference-beam autocollimator with nanoradian sensitivity from mhz to khz and dynamic range of 10(7.). *The Review of scientific instruments*, 84:095007, 09 2013.